


- Manufactured in a certified ISO Class 7 cleanroom.
- Made entirely of 100% USP Class VI components from the United States.
- Assemblies do not use bonding agents and are designed for one continuous flow path.
- Gamma Irradiated in the United States.

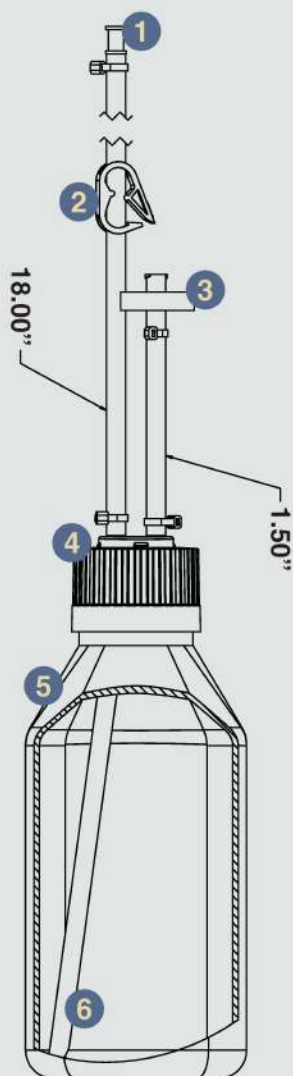

Foxx Life Sciences' line of EZbio® Single Use Sterile Media Bottles with Tubing and Vent Filters are perfect for fluid transfer and storage. Assemblies are made entirely with USP Class VI components and Gamma Irradiated in the United States. Assemblies take advantage of Foxx's patented VersaCap® technology which means there is no machining necessary to produce the product. Perfect for single use applications and to phase out costly stainless steel apparatuses that are expensive to clean and validate.

Foxx also offers the ability to provide and produce custom components, assemblies, and variations in a timely manner. For more information please contact us at sales@foxxlifesciences.com or by phone at 1-603-890-3699.

- 1 Hose barb plug keeps line closed until it is removed for use.
- 2 Pinch clamp that easily allows the user to control flow into the container, or to close off line entirely.
- 3 0.2µm Hydrophobic PTFE filter allows for ventilation while preventing media from escaping in the event of a spill.
- 4 VersaCap® technology ensures a closed system without any machining residue or bonding agents.
- 5 Media bottles are available in five different sizes: 125ml, 250ml, 500ml, 1L, 2L.
- 6 Rigid dip tube with notch to allow for easy filling and extraction of media.

**Available in both
PC and BPA-Free
PETG!**

**MADE
IN** 

ISO 9001 • ISO 13485





	DESCRIPTION	VOLUME	PART NUMBER
PETG	EZbio Single Use Kit, Vented VersaCap 38-430	125mL	437-0003-FLS
	EZbio Single Use Kit, Vented VersaCap 38-430	250mL	437-0004-FLS
	EZbio Single Use Kit, Vented VersaCap 38-430	500mL	437-0005-FLS
	EZbio Single Use Kit, Vented VersaCap 38-430	1000mL	437-0006-FLS
	EZbio Single Use Kit, Vented VersaCap 53B	2000mL	437-0007-FLS

	DESCRIPTION	VOLUME	PART NUMBER
PC	EZbio Single Use Kit, Vented VersaCap 38-430	125mL	438-0003-FLS
	EZbio Single Use Kit, Vented VersaCap 38-430	250mL	438-0004-FLS
	EZbio Single Use Kit, Vented VersaCap 38-430	500mL	438-0005-FLS
	EZbio Single Use Kit, Vented VersaCap 38-430	1000mL	438-0006-FLS
	EZbio Single Use Kit, Vented VersaCap 53B	2000mL	438-0007-FLS



ISO 7 Cleanroom

Foxx Life Sciences has a certified ISO 7 Cleanroom with 6 full assembly lines just 32 miles north of Boston, MA! We provide low cost, rapid assemblies for custom projects while maintaining stringent quality standards. All assemblies are built following our ISO 9001 and 13485 quality procedures. In addition, Foxx offers package validation for new products along with complete material sourcing. Contact us today to find out more about all of our capabilities and opportunities that we can offer you!