



VELTEK ASSOCIATES, INC.

TECHNICAL DATA FILES



Core2Write® Labels Blank & Custom Cleanroom Labels

10APR2024

Copyright © 2024 Veltek Associates, Inc.

All Rights Reserved. Confidential and Proprietary Property of Veltek Associates, Inc. No part of this document may be reproduced or disclosed, in whole or in part, without the specific written permission of Veltek Associates, Inc.



Core2Write® Labels

Blank & Custom Cleanroom Labels

PRODUCT OVERVIEW

Core2Write cleanroom labels are made from VAI's cellulose free CleanPrint 10[®] synthetic writing substrate. Due to CleanPrint 10's cellulose free nature, using **Core2Write** labels in the core and critical environments considerably reduces fibers and particulates and potential sources of cross-contamination. Blank or custom labels come standard with cleanroom compatible adhesives and cores. CleanPrint 10 is smear resistant so it is compatible with chemical disinfecting, is autoclavable, and is exceptionally durable. Labels are available sterile and non-sterile, on rolls or sheets, and in numerous sizes and color combinations.

In addition to custom labels, VAI[®] offers standard blank white **Core2Write** labels that are ideal for facility in-house printing. These labels are either available on sheets or in rolls, as sterile or non-sterile, and have standard bagging and case quantities.

Sterile **Core2Write** labels are gamma irradiated to 10⁻⁶ sterility assurance level and have been completely validated for sterility and a ten-year shelf life. Lots are processed in a clean environment, completely traceable, and are sterility tested according to current USP compendium. Sterile **Core2Write** labels are delivered each time with lot specific documentation including a Certificate of Sterility and Certificate of Irradiation.

Core2Write custom labels are specifically made for your organization's needs. Depending on requirements, the following items can be customized: color, size, company name, logo, SOP name and number, document revision number, custom/unique numbering, RFID integration, bar coding, QR coding, temperature sensitivity, specific entry requirements, bagging, transportation requirements, labeling, and sterility. Custom validation packages are also available.



Features and Benefits

- Fully customizable – numerous sizes and colors available
- Blank labels available for in-house facility printing
- Rolls or sheets available
- Printed on cellulose free CleanPrint 10
- Reduces fibers, particulates, and potential contaminants in the core and critical environments
- Smear resistant, chemical compatible, and exceptionally durable
- Standard 7 mil thickness
- Cleanroom compatible adhesive and core
- Sterile labels are gamma irradiated at 10⁻⁶ SAL
- Validated for sterility and shelf life with either standard or custom validations available
- 10-year expiration dating
- Completely lot traceable with lot specific documentation
- Sterile products are lot sterility tested according to current USP compendium

Patents: www.sterile.com/patents




Core2Write® Labels

Product Uses & Ordering Information

Product Uses

Core2Write cleanroom labels can be used to:

- Document GMP operations
- Label or tag GMP containers
- Label or tag transfer cans/tanks/bottles/bags
- Denote clean or needs cleaning
- Denote waste, experimental, reject, quarantine, in-process, and many more
- Use in grade A-D environments
- With RFID tracking systems

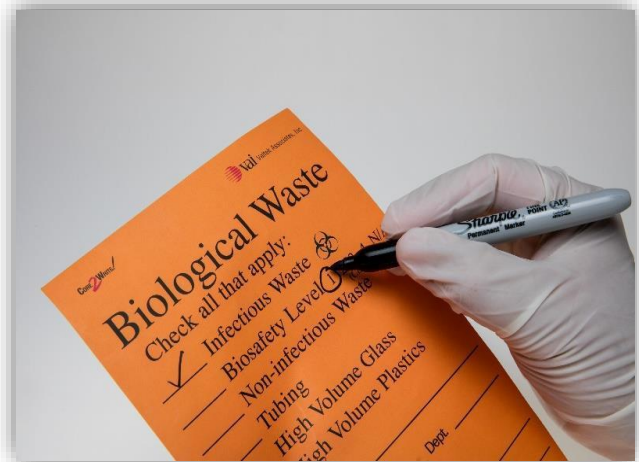

Core2Write labels can have RFID inlays embedded into the material so that assets, equipment, and procedures can be tracked and found with a RFID reader. Core2Write products, via their RFID inlay capabilities, are compatible with VAI's Core2Scan end-to-end tracking system. Core2Scan is a total asset tracking and asset management solution that incorporates hardware, proprietary software, and RFID compatible labels into a complete package to alleviate tracking and traceability challenges.

ORDERING INFORMATION

Custom Labels

The Core2Write label product line starts with a customer evaluation of which labels are required and to what specifications. Once determined, labels and artwork are digitally designed into the product required and printed. Custom VAI part numbers that are specific to your firm will be created for re-ordering purposes.

*Standard Blank Labels**

Part Number	Description	Qty/Case
C2WR-2X1-BLK-01	Core2Write Label, Rolls, 2x1, Blank, Non-Sterile	1000 labels/roll, 4 rolls/case, 4000
C2WR-2X1-BLK-02	Core2Write Label, Rolls, 2x1, Blank, Sterile	1000 labels/roll, 4 rolls/case, 4000
C2WR-3X1-BLK-01	Core2Write Label, Rolls, 3x1, Blank, Non-Sterile	1000 labels/roll, 4 rolls/case, 4000
C2WR-3X1-BLK-02	Core2Write Label, Rolls, 3x1, Blank, Sterile	1000 labels/roll, 4 rolls/case, 4000
C2WR-3X1.5-BLK-01	Core2Write Label, Rolls, 3x1.5, Blank, Non-Sterile	1000 labels/roll, 4 rolls/case, 4000
C2WR-3X1.5-BLK-02	Core2Write Label, Rolls, 3x1.5, Blank, Sterile	1000 labels/roll, 4 rolls/case, 4000
C2WR-4X2-BLK-01	Core2Write Label, Rolls, 4x2, Blank, Non-Sterile	1000 labels/roll, 4 rolls/case, 4000
C2WR-4X2-BLK-02	Core2Write Label, Rolls, 4x2, Blank, Sterile	1000 labels/roll, 4 rolls/case, 4000
C2WR-4X3-BLK-01	Core2Write Label, Rolls, 4x3, Blank, Non-Sterile	1000 labels/roll, 4 rolls/case, 4000
C2WR-4X3-BLK-02	Core2Write Label, Rolls, 4x3, Blank, Sterile	1000 labels/roll, 4 rolls/case, 4000
C2WR-2X1-SHEET-01	Core2Write Label, Sheet, 2x1, Blank, Non-Sterile	500 sheets/case, 36 labels/sheet, 18000
C2WR-2X1-SHEET-02	Core2Write Label, Sheet, 2x1, Blank, Sterile	500 sheets/case, 36 labels/sheet, 18000
C2WR-3X1-SHEET-02	Core2Write Label, Sheet, 3x1, Blank, Sterile	500 sheets/case, 20 labels/sheet, 10000
C2WR-4X2-SHEET-02	Core2Write Label, Sheet, 4x2, Blank, Sterile	500 sheets/case, 10 labels/sheet, 5000
C2WR-4X3-SHEET-02	Core2Write Label, Sheet, 4x3, Blank, Sterile	500 sheets/case, 6 labels/sheet, 3000

*Standard blank labels are also available in non-sterile, sheet label options, and in additional sizes. Contact your local sales representative for more information.